OxyGold™ G Operating Instructions

Introduction

These operating instructions refer to Hamilton's OxyGold G.

Name	Ref
OxyGold G 120	237395
OxyGold G 225	237396

Hamilton sensors are quality products manufactured with the latest technology. You can only attain the highest accuracy and a maximum durability by exactly following the instructions below.

These instructions should be read, understood and followed by all staff using the device. Hamilton can assume no responsibility for damage and operational disruptions arising from failure to observe these instructions.

Liability

The liability of Hamilton Bonaduz AG is detailed in the document «General Terms and Conditions of Sales and Delivery (GTS)», chapter 12.

Hamilton expressly shall not be liable for direct or indirect losses arising out of the utilization of the sensors. It must in particular be ensured in this conjunction that malfunctions can occur on account of the inherently limited useful life of sensors contingent upon their relevant applications. The user is responsible for the calibration, maintenance and punctual replacement of the sensors. In the case of critical sensor applications, Hamilton recommends using redundant measurement points in order to avoid consequential damages. The user shall be responsible for taking suitable precautions in the event of a sensor failure.

Intended use

The OxyGold electrode is designed to measure the partial pressure of dissolved oxygen in liquids. The sensors is optimized for the measurement of trace quantities of dissolved oxygen. Its lower detection limit is 1 ppb. Typical application is the measurement of boiler feed water.

The OxyGold G sensor has a built- in temperature sensor (NTC 22 kOhm). This temperature sensor is to be used only to compensate the oxygen signal and not to control process temperature.

If the sensor is used in potentially explosive atmospheres, refer to the sections «Instructions for use in potentially explosive atmospheres».


⚠ Safety aspects

This sensor is to be used only as intended and is to be kept in a condition that ensures complete safety. The specifications given in the section «Technical Data» as regards temperature, pressure etc. may under no circumstances be exceeded. Inappropriate use or misuse can be dangerous. Assembly and maintenance may only be done by trained personnel.

Since the inner part of the sensor is glass, handle it with care. The top of the sensor can be damaged by impacts. When changing the electrolyte it is recommended to wear safety glasses and protective gloves. Contact of the electrolyte with skin, eyes or mucous membrane is to be avoided (see Safety Data Sheet Ref 608929).

If mechanical damage to the sensor occurs, electrolyte may leak out. Ensure that when screwed into the process the PG13,5 thread and the O-ring incur no damage. O-rings are consumable items and must be replaced regularly – at least once a year. Even where all necessary safety measures have been complied with, there is still a further risk involving leaks or mechanical damages, where there are seals or screws, gases or liquids could leak out undetected. Before dismantling the sensor, make sure that there is no over-pressure in the process.

Every product for shipment or sent back for repair must be decontaminated. If working with hazardous liquids observe and carry out the maintenance procedures, paying particular attention to cleaning and decontamination. If the product becomes contaminated with biohazardous, radioactive or chemical material, it should be cleaned.

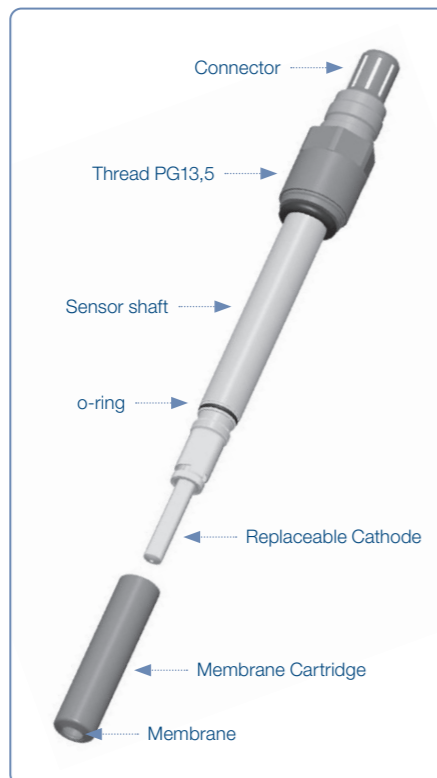First inspection

During unpacking please check for possible mechanical defects. If there is any problem please return the sensor in its original packing to Hamilton quoting the invoice number and serial number (on the PG13,5 thread part).

Initial Operation

- The sensor is packed up dry. Remove the membrane cartridge rinse it once with the electrolyte Oxylyte G. Then pour in 1.5 ml of electrolyte using the enclosed plastic pipette. Carefully screw the membrane cartridge back in.

⚠ **Attention**
Do not touch or bend the glass body. The electrolyte is aggressive.
- Screw the sensor into the housing with the PG thread. Ensure that the torque of 1.5 Nm will not be exceeded to avoid damaging of the O-ring. The membrane must be pointing downwards.
- Connect the OxyGold with the amplifier. Its polarization voltage should be -670 mV. After a minimum period of 2 hours, the sensor is ready for calibration. This stabilization time is necessary to reach the specification of the zero current.
- To calibrate, follow the instrument manual. Air/gas bubbles may stick to the sensitive area of the sensor. As a consequence, the measurement value might be wrong. In case of a shutdown of the power supply (230 V) the measurement value could be wrong.

Dismounting the sensor

Before removing the sensor from the measuring setup, always make sure that the setup is pressureless and cold and that no process medium can be accidentally spilled. Unscrew the PG13,5 thread and pull out the sensor. Only turn the PG13,5 thread, not the sensor body, otherwise you unscrew the cap.

If working with hazardous liquids observe and carry out the maintenance procedures, paying particular attention to cleaning and decontamination.

Cleaning, assembly and maintenance should be performed by personnel trained in such work. Do not use any abrasive tissues or cleaning materials and do not use any cleaning chemicals other than described above.

When removing and cleaning the sensor, it is recommended to wear safety glasses and protective gloves.


VP Connector

When the OxyGold is being disassembled from the cable, there is some risk of water entering into the VP female on the cable side. Ensure that the cable/connector is dry when unscrewing the electrode from the cable.

Sensor check and maintenance

For a periodical check of the OxyGold sensor please proceed as follows:

- Hold a stabilized sensor in air and read measurement value. Hold the sensor in nitrogen flow and read measurement value (e.g. place nitrogen in a plastic bag). Measurement value must be under 2% of air value. Too high a current in nitrogen is due to an exhausted electrolyte, or a defective membrane. In such a case replace the membrane cartridge as well as the electrolyte solution.
- If there is a sluggish response and/or too small a current in air, follow these steps:

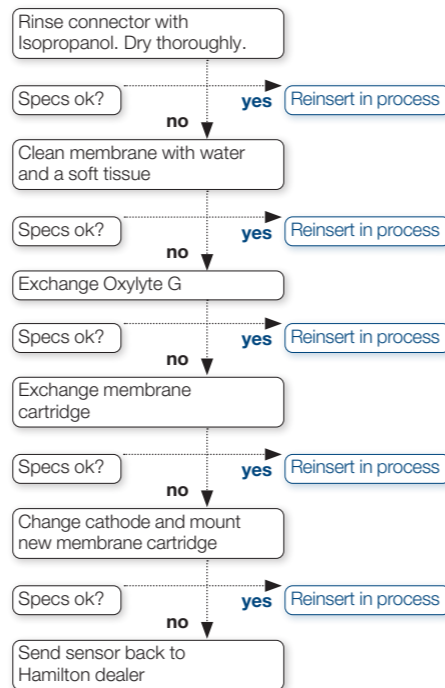«Specs OK» means: current in air of stabilized sensor at 25°C is 200 – 500 nA. $t_{98\%} < 60$ s at 25°C.

Replacing electrolyte and membrane cartridge

The OxyGold membrane cartridge is replaced as follows:


- Put the sensor in an upright position and carefully unscrew the membrane cartridge (E), empty and rinse with water.

⚠ **Attention**
Do not touch or bend the glass body. The electrolyte is aggressive.
- Check the small O-ring (D) above the cathode. Replace it, if it is damaged.
- Use the 1 ml plastic pipette of the Membrane Kit to dispense 1.5 ml of electrolyte Oxylyte G (Ref 237139) into the new membrane cartridge.
- Carefully screw the cartridge onto the sensor shaft. Any spillage of electrolyte should be rinsed away with water.
- Unscrew and empty membrane cartridge. Repeat steps 3 and 4.

Replacing the Cathode

If the cathode is changed, the membrane must also be changed to get the full performance of the sensor.

- Put the sensor in an upright position and carefully unscrew the membrane cartridge (E).
- Flush the replaceable cathode with deionized water then dry the metal parts. Do not touch anode and cathode.
- Hold the cathode on the metal part (B) in front of the thread (A) and the sensor on the shaft (C) and pull apart. Just pull, do not turn!
- Please check that all contacts are dry and clean.
- Replace cathode for a new (Ref 237306 (blue) for OxyFerm; Ref 237427 (yellow) for OxyGold B; Ref 237437 (black) for OxyGold G). Rotate the cathode until the correct position is found and then push the cathode into the shaft.
- Follow now the instructions of «Replacing electrolyte and membrane cartridge».

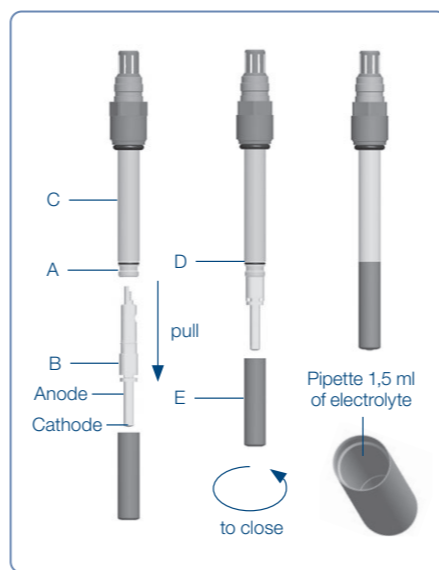
Electrical connections

Please use the appropriate Hamilton VP-cables. See Ordering Information.

Sensor connector:

A	cathode
B	anode
C	-
D	contact to metal shaft*
E, F	temperature

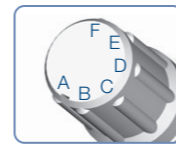*To improve the signal stability the pin D of the VP connector can be grounded separately.

ATTENTION: This grounding must be galvanically isolated from the anode and cathode, otherwise the sensor lifetime will be shortened significantly!

Disposal

The design of Hamilton sensors minimizes environmental impact. According to the EU directive 2012/19/EU the Hamilton sensors should be disposed of waste of electrical and electronic equipment, and not in municipal waste. They can also be sent back to Hamilton for disposal.

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Instructions for use in potentially explosive atmospheres

ATEX / UKEX / IECEx marking:

Gas: CE 0035 / UKEX 2571
II 1/2 G Ex ia IIC T4/T5/T6 Ga/Gb

Dust: CE 0035 / UKEX 2571
II 1/2 D Ex ia IIC T x °C Da/Db

Manufacturer:
Hamilton Bonaduz AG, 7402 Bonaduz, Switzerland

EU-Type Examination Certificate: TÜV 03 ATEX 7005 X
UK Type Examination Certificate: TÜV 21 UKEX 7047 X
IECEx Certificate of Conformity: IECEx TUR 14.0001 X

EU-Type Examination Certificate, UK Type Examination Certificate, and IECEx Certificate of Conformity can be downloaded from www.hamiltoncompany.com.

Conditions described in the ATEX EU-Type Examination Certificate, UK Type Examination Certificate, or IECEx Certificate of Conformity must be adhered to.

⚠ Attention

In case a gas atmosphere and a dust atmosphere are or could be present at the same time, the risk of explosion must be examined carefully and special precautions may be necessary.

Coupling situation

- The operator of equipment in potentially explosive atmospheres is responsible for ensuring that all components of the system are certified for that area classification and are compatible with each other.
- In gas atmospheres, the Ex approval is not dependent on any conditions regarding mounting of the sensor. Use with particular Hamilton Armatures can however lead to an increase in the allowable power of the transmitter. See section «Connection to

Transmitter». In dust atmospheres, special restrictions may apply when using armatures made of plastic material.

- O-rings having a sealing function between Ex-Zone 0 and 1, 20 and 21 respectively must be replaced at each dismantling of the sensor.
- Sensors, transmitter and other required equipment are to be set up within a single potential equalizing system.
- When the sensor is assembled, the attached ATEX / UKEX / IECEx sticker is to be attached to the cable in an easily visible place, as near as possible to the sensor itself (this sticker is to indicate that an ATEX / UKEX / IECEx approved sensor is in use). The sticker should not be removed at any time.

⚠ Attention

Inappropriate fitting or dismantling of the armature or the sensor may lead to unintended release of an explosive atmosphere.

Connection to Transmitter

OxyGold G is suitable for connection to an intrinsically safe electrical circuit with protection level ia. The operator of the equipment must ensure that the allowable electrical values for the sensor all exceed those for the transmitter. None of the values cited for electrical power, voltage and current may be exceeded in total (measurement and temperature circuit together). The electrical values for OxyGold G are as follows:

U =24 V; I = 173 mA; P: Function of process temperature.

If the OxyGold G sensor is used with a RetractoFit or FlexiFit armature as per the table below, the values of type 5 are the correct ones. In all other coupling situations the values of type 4 are the correct ones.

Coupling situation for type 5	Sensor a-length 120	Sensor a-length 225
RetractoFit Serie*	–	yes
RetractoFit Bio Serie*	–	yes
FlexiFit Serie*	yes	–
Retractex Serie*	–	yes

* Several Ref possible.

For gas atmospheres:

The maximum allowable electric power P of the transmitter (measurement and temperature circuit together) depends on the desired process temperature and the coupling situation:

Type 5	P ≤ 250 mW
T4	125 °C
T5	90 °C
T6	75 °C

Type 4	P ≤ 30 mW	P ≤ 60 mW	P ≤ 90 mW
T4	129 °C	129 °C	128 °C
T5	94 °C	94 °C	93 °C
T6	68 °C	57 °C	46 °C

For dust atmospheres:

Table for the calculation of the maximum surface temperature «x» of the sensor as a function of the maximum electrical power of the transmitter «P» and ambient/process temperatures «T_a». «x» must be smaller than the ignition temperature of dust involved.


Type 5	P ≤ 250 mW
	x=T _a +5 °C

Type 4	P ≤ 30 mW	P ≤ 60 mW	P ≤ 90 mW
	x=T _a +12 °C	x=T _a +23 °C	x=T _a +34 °C

The temperature limit (130°C) must always be respected.

Technical Data

Measurement Principle	Electrochemical reduction of oxygen
Measuring Range	1 ppb to 40 ppm (DO) or 0.02 to 1000 mbar (pO ₂)
Oxygen Consumption	ca. 100 ng/h in air at 25 °C
Required Flow	≥ 0.1 m/s
Response Time t _{98%}	30 to 60 s at 25 °C, from air to nitrogen
Temperature Sensor	NTC 22 kOhm
Current in air	see fig.

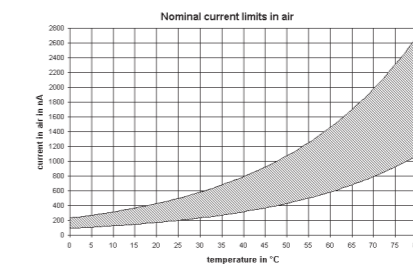Drift at Room Temp. under Constant Conditions	< 1 % per week
Pressure Range	0 to 12 bar
Max. CO ₂ partial pressure	0.01 bar
Residual current of current air	< 0.02 %
Membrane / Cap	Delivered with Optiflow
Shaft diameter	12 mm
Mounting	PG13,5
Wetted Parts	Stainless steel 1.4435, Silicone, EPDM
Polarization Voltage and Time	-670 ± 50 mV, ≥ 2 hours
Operating Temperature Range	0 to 130 °C
Temperature Response	3.1 %/K
Storage Temperature	-10 to 60 °C

Ordering information

OxyGold G 120	Ref 237395
OxyGold G 225	Ref 237396

Accessories

VP-Cable, 1 m	Ref 355108
VP-Cable, 3 m	Ref 355109
VP-Cable, 5 m	Ref 355110
VP-Cable, 10 m	Ref 355111
VP-Cable, 20 m	Ref 355112
OxyGold Membrane Kit	Ref 237135
Replacement Cathode OxyGold G	Ref 237427
Oxylyte G Electrolyte 50 ml	Ref 237139
Polarization module G for OxyGold G	Ref 237350

Important notice

The contents of these operating instructions are subject to change without notice. Technical changes reserved. The above notwithstanding, Hamilton Bonaduz AG can assume no responsibility for any errors in these operating instructions or their consequences.

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© 2023 Hamilton Bonaduz AG. All rights reserved.